


SUCCESSFUL WOUND TREATMENT STARTS WITH YOU

Your healthcare provider has recommended Reaplix 3C Patch® treatment for you.

3C Patch® draws on your body's capacity to heal. In short, you aid your own healing process.

A highly effective treatment, 3C Patch® is proven to substantially increase your chances of complete wound closure.

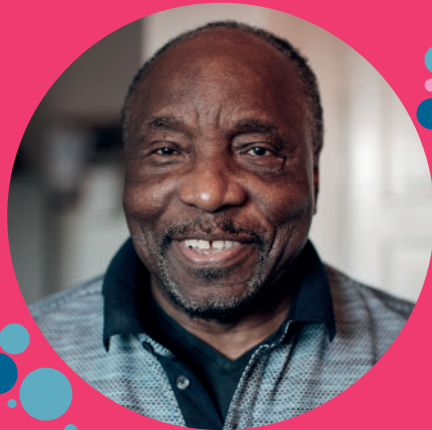

IMPROVED CHANCES OF HEALING¹

With 3C Patch®, your wound is 89% more likely to heal compared with best standard of care.

89%

1. Game F et al, LeucoPatch system for the management of hard-to-heal diabetic foot ulcers in the UK, Denmark, and Sweden: an observer-masked, randomized controlled trial. The Lancet Diabetes Endocrinology 2018



OUR STORY

Reaplix has established that successful wound treatment can be accomplished by using the body's own capacity to heal.

3C Patch® is a highly effective wound therapy that draws on your body's unique ability to aid its own healing process, and is recommended by the International Working Group on Diabetic Foot (IWGDF) Guidelines.

CONTACT

USA

REAPPLIX INC / US
info@reaplix.com
+1 866 327 2824

EUROPE (HQ)

REAPPLIX A/S HQ / EUROPE
info@reaplix.com
+45 88 16 80 66

3C Patch® is formerly known as LeucoPatch. 3C Patch is a trademark of Reaplix A/S, registered in the EU and the US. 3CP is a trademark of Reaplix A/S, registered in the EU and pending in the US. 00950 EU&US ver. 4.0



EVERY WOUND IS PERSONAL

YOUR 3C PATCH® TREATMENT EXPLAINED



YOUR BODY HOLDS THE KEY TO HEALING


The 3C Patch® is made entirely from 18 ml of your own blood – nothing else. Your blood sample is placed in an automated centrifuge, where it's processed to generate a three layer patch.

Each 3C Patch® takes approximately 15-20 minutes to produce. If you are taking blood-thinning medication, the process time may be longer.

The 3C Patch® is applied directly to your cleaned wound where it remains until your next treatment the following week. If your wound is larger than 5 cm², you may require more than one patch.



YOUR 3C PATCH® HAS THREE LAYERS

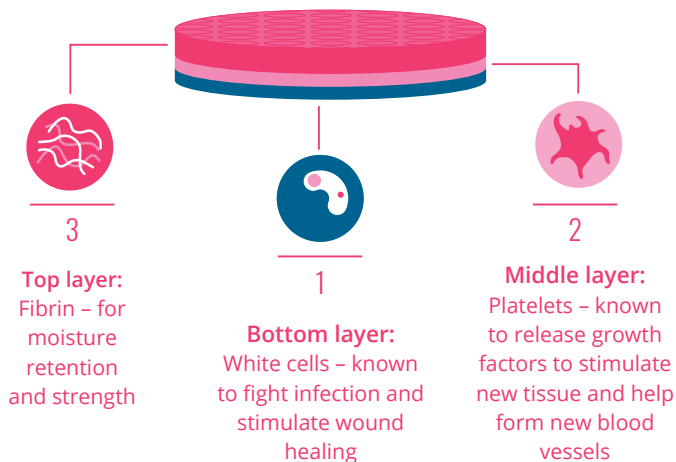

WHAT YOU WILL EXPERIENCE WITH THE 3C PATCH® TREATMENT

The 3C Patch® usually dissolves in a week. If there is any loose material in the wound at your next treatment, your healthcare provider will remove it before applying a new patch.

The number of treatments you require will depend on the size and depth of your wound and how you respond to treatment.

After two to three treatments, you should see an effect. Eight to ten treatments are not uncommon.

HOW TO CARE FOR YOUR 3C PATCH®

A primary dressing is applied to keep the 3C Patch® in place. It must not be removed until your next 3C Patch® treatment. Please share this information with those helping you with dressing changes.

A secondary dressing on top of the primary dressing can be changed as often as required. When showering or bathing, please protect the dressing with a plastic bag to keep it dry.

It's important that you follow your healthcare provider's recommendations for wound care and supportive footwear.

WHAT TO LOOK OUT FOR

Any open wound is at risk of infection. During the course of treatment, it's important to be aware of any changes to the wound, which may include:

- REDNESS
- SWELLING
- HEAT
- PAIN

IF YOU EXPERIENCE ANY OF THESE SYMPTOMS, PLEASE CONTACT YOUR HEALTHCARE PROVIDER AS SOON AS POSSIBLE.